

FROM ALLESEE ORTHODONTIC APPLIANCES

Active Plates

Schwartz and Sagittal Appliances (Active Plates) can be custom designed for many different functions and combinations of functions. Screws can be placed in transverse and Sagittal positions or used for single tooth movement. Special screw types are available for specific areas of the arch such as just the anterior and not the posterior. Posterior acrylic coverage is available as a means to disocclude the opposing arch and aid in the movement of these covered teeth.

Other popular options include finger springs (either acrylic borne or soldered), habit cribs, pearls, headgear tubes and elastic hooks.

Additional information is available upon request to our Customer Service number:

1-800-262-5221

Fax 1-414-886-6879


www.ormco.com/aoaapps.html




P.O. Box 725
Sturtevant, WI 53177

AOA is a subsidiary of Ormco Corporation


Rev. B

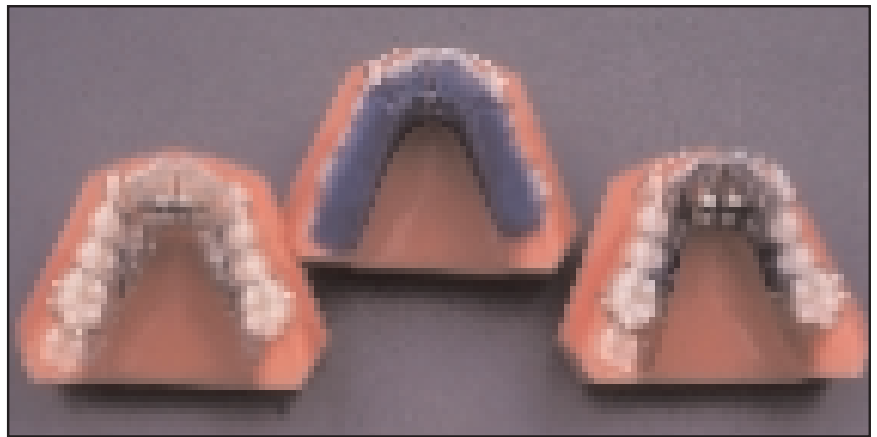


Lower Schwartz Appliances with and without labial bow.

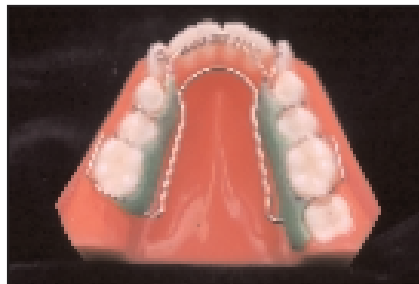


Lower Jackson Appliance for a persistent, gentle expansion force.

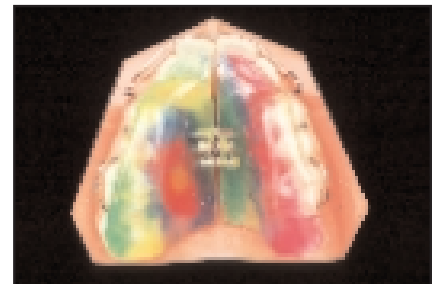

Upper Schwartz Appliance with occlusal coverage and AOA's standard labial bow and clasp.

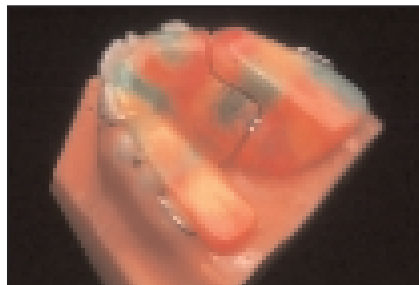

Upper and lower Sagittal Appliances. Occlusal coverage can be added upon request.

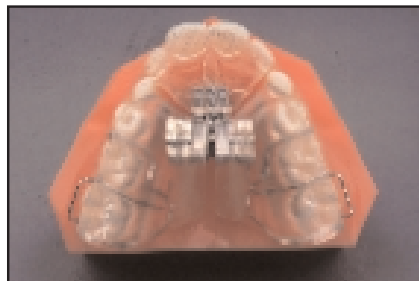

3-Way Expanders. Available with three independent jackscrews or a single screw with both transverse and A-P expansion capability.



Active plate with either stainless steel or TMA rectangular distalizing springs. Short anterior bite plane to disocclude the posterior teeth.